

TEOS characterization of 2D materials – from graphene to TMDCs

Marc Chaigneau¹

Yoshito Okuno¹, Andrey Krayev², Filippo Fabbri³

¹HORIBA Scientific, passage Jobin Yvon, Palaiseau, France

²AIST-NT Inc., Bel Marin Keys Blvd, Novato, USA

³MEM-CNR Institute, Parco Area delle Scienze, Italy

marc.chaigneau@horiba.com

Raman spectroscopy has become a leading technique in graphene and 2D materials research. It is widely used to determine the number and orientation of layers, the quality of the materials, the presence of defects and their density, and the effects of strain or doping elements. Tip Enhanced Raman Scattering (TERS) is today a mature technique that gives access to all those mentioned information down to the nanoscale, thus providing important insights into the structure of such 2D materials.

In this talk, we'll report results of TEOS (Tip Enhanced Optical Spectroscopy, including TERS and TEPL, Tip Enhanced Photoluminescence) characterization of graphene, functionalized graphene oxide and 2D semiconductors, MoS₂ and WS₂. We will report on how TEOS characterization can reveal very localized spectra change correlated with the number of layers, the local position of the Fermi level, and breaks of translation symmetry.

We will discuss the impact of TEOS characterization in the studies of single-layer and few layer 2D materials and emphasize the capability of seeing deep-sub-diffraction-limited details that can be today correlated with complementary Kelvin probe microscopy and photo-current techniques.

Figures

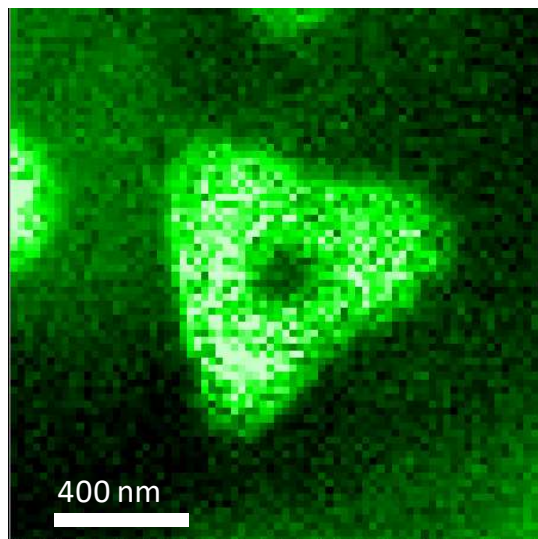


Figure 1: TERS/TEPL image – MoS₂ flake on Si

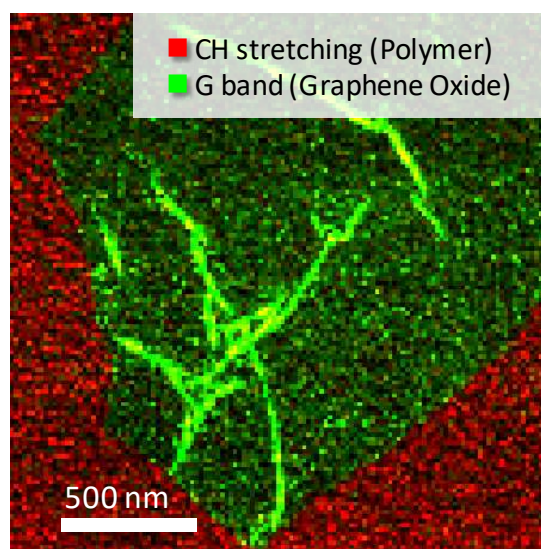


Figure 2: TERS image – Graphene oxide on gold

